



Safety & user manual

4-POINT HIP BELT

Articles: 6.70042.X, 6.70043.X,
6.70044.X, 6.70045.X, 6.70047.X.

KBSUM-4HB0925



This manual provides essential information for the safe use and maintenance of the 4-point hip belt range manufactured by Kinetic Balance. It is important that both the user and/or caregiver read and fully understand these instructions before use.

Intended use:

⚠ This product is exclusively designed for positioning the pelvis/hip of individuals seated in wheelchairs or other seating systems. It must **not** be used as a transportation belt or as a restraint device. Improper use may result in serious safety risks.

⚠ The belt should always be securely fitted around the pelvis or over the thighs during use. Installation and adjustments must be carried out by a qualified seating specialist.

⚠ If a user's physical or cognitive abilities limit their ability to operate the belt safely, a caregiver must be present at all times. Caregivers should be properly instructed in fastening and releasing the product. This product does not replace the need for caregiver supervision.



Periodic Safety and Maintenance checks:

To ensure user safety, regularly inspect the following components:

- Proper functioning of the buckle
- Signs of wear or damage to padding, stitching, or webbing
- Correct positioning of the product on the user

If any malfunctioning or wear is found on the belt, please contact your supplier for repair or replacement.

▲ In the event of a serious incident, report it to the manufacturer and the relevant local authority.

Cleaning instructions:

Hand wash or machine wash at 30°C. Do not bleach, tumble dry, or iron. To protect the product during washing, use a laundry bag. If required for more intensive cleaning, the product may be machine washed at 60°C.



Recycling:

At the end of its lifecycle, the product must be dismantled and disposed of in accordance with applicable local regulations for each material.

Warranty:

The product includes a 12-month warranty covering manufacturing defects.

Mounting Kit:

SIZE	WEBBING CLIP			
	15mm	20mm	25mm	38mm
XS		2pc	2pc	
S			2pc	2pc
M			2pc	2pc
L			2pc	2pc
XL			2pc	2pc

SIZE	CAM BUCKLE			
	15mm	20mm	25mm	38mm
XS			4pc	
S			2pc	2pc
M			2pc	2pc
L			2pc	2pc
XXL			2pc	2pc



SIZE	SLIDER BUCKLE			
	15mm	20mm	25mm	38mm
XS		2pc	2pc	
S			2pc	2pc
M			2pc	2pc
L			2pc	2pc
XXL			2pc	2pc




Mounting a 4-point hip belt

The mounting of a 4-point hip belt should be determined based on the clinical goal. Always consult a seating specialist to determine optimal positioning and belt angle.

Installation steps:

1. Identify the appropriate mounting points for the primary and secondary straps on the wheelchair or seating system
2. Attach the mounting hardware securely
3. Center the belt on the user and thread the webbing through the mounting
4. Adjust the belt length
5. Tighten the belt to ensure a snug and secure fit

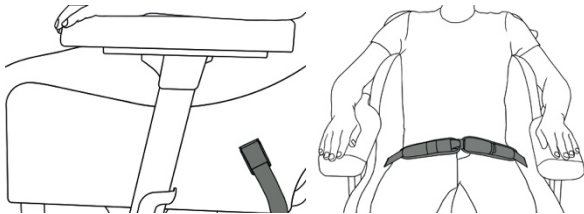

Want to know more about positioning?
Watch our e-learnings.

E-learnings



Adjustment of the belt


During use the belt can be adjusted using the D-ring(s). Ensure the buckle remains centered on the user during adjustments.

Adjustments should only be made to:

- Facilitate easy fastening and unfastening
- Ensure the belt remains snug



Fastening and Loosening buckles

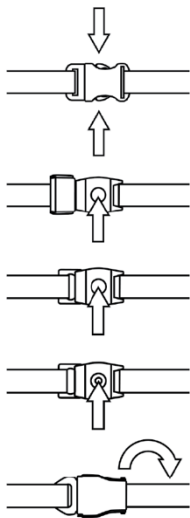

Symbols list



The product is manufactured by:
Kinetic Balance



Address:
Scheepmakerstraat 27
2222AB
Katwijk aan Zee
The Netherlands

Email: mail@positioningproducts.nl
Phone: +31 (0)71-2048424
Website: www.positioningproducts.nl
Chambre of commerce: 27339364

