



Safety & user manual

ANKLE POSITIONING

Articles: 8.70001.X, 8.70002.X

KBSUM-ANP0925



This manual provides essential information for the safe use and maintenance of the ankle positioning range manufactured by Kinetic Balance. It is important that both the user and/or caregiver read and fully understand these instructions before use.

Intended use:

⚠ This product must only be used in combination with a hip belt.

⚠ This product is exclusively designed for positioning the feet on the footrest of individuals seated in wheelchairs or other seating systems. It must not be used as a transportation belt or as a restraint device. Improper use may result in serious safety risks.

⚠ The product should always be securely fitted around the ankles during use. Installation and adjustments must be carried out by a qualified seating specialist.

⚠ If a user's physical or cognitive abilities limit their ability to operate the belt safely, a caregiver must be present at all times. Caregivers should be properly instructed in fastening and releasing the product. This product does not replace the need for caregiver supervision.



Periodic Safety and Maintenance checks:

To ensure user safety, regularly inspect the following components:

- Proper functioning of the buckles or hook and loop closure
- Signs of wear or damage to padding, stitching, or webbing
- Correct positioning of the product on the user

If any malfunctioning or wear is found on the belt, please contact your supplier for repair or replacement.

▲ In the event of a serious incident, report it to the manufacturer and the relevant local authority.

Cleaning instructions:

Hand wash or machine wash at 30°C. Do not bleach, tumble dry, or iron. To protect the product during washing, use a laundry bag. If required for more intensive cleaning, the product may be machine washed at 60°C.

Recycling:

At the end of its lifecycle, the product must be dismantled and disposed of in accordance with applicable local regulations for each material.



Warranty:

The product includes a 12-month warranty covering manufacturing defects.

Mounting Kit:

SIZE	METAL HANDLE
XS	4pc
S	4pc
M	4pc
L	4pc

SIZE	SLIDER BUCKLE			
	15mm	20mm	25mm	38mm
XS		2pc		
S			2pc	
M				2pc
L				2pc




Mounting an Ankle Positioning

The mounting of the ankle positioning should be determined based on the clinical goal. Always consult a seating specialist to determine optimal positioning and angle.

Installation steps:

1. Locate the mounting points on the footplate
2. Attach the metal brackets securely to these points
3. Thread the webbing through the slider buckle, then loop it around the metal bracket and pass it back through the slider buckle
4. Position the ankle support around the ankle, ensuring it sits over the malleoli. The buckles may be positioned on either the inside or outside of each foot

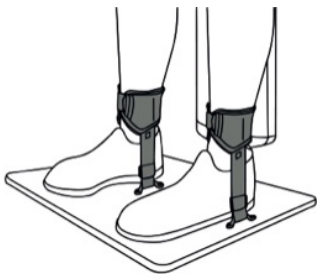

Want to know more about positioning?
Watch our e-learnings.



E-learnings



Adjustment of the abduction positioning

Daily adjustments may be needed because of variations in footwear such as socks, shoes and slippers.

- For hook-and-loop closures, daily adjustment is achieved by tightening the closure until a secure and snug fit is obtained
- For side-release buckles, use the webbing adjustment on the buckle to accommodate daily variations

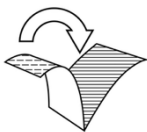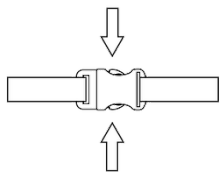

Symbols list



The product is manufactured by:
Kinetic Balance



Address:
Scheepmakerstraat 27
2222AB
Katwijk aan Zee
The Netherlands

Email: mail@positioningproducts.nl

Phone: +31 (0)71-2048424

Website: www.positioningproducts.nl

Chambre of commerce: 27339364

