



Safety & user manual

CHEST BELT

Articles: 400X, 402X

KBSUM-CB0925



This manual provides essential information for the safe use and maintenance of the chest belt range manufactured by Kinetic Balance. It is important that both the user and/or caregiver read and fully understand these instructions before use.

Intended use:

⚠ This product is exclusively designed for positioning the trunk of individuals seated in wheelchairs or other seating systems. It must **not** be used as a transportation belt or as a restraint device. Improper use may result in serious safety risks.

⚠ This product should always be securely fitted around the trunk during use. Installation and adjustments must be carried out by a qualified seating specialist.

⚠ If a user's physical or cognitive abilities limit their ability to operate the belt safely, a caregiver must be present at all times. Caregivers should be properly instructed in fastening and releasing the product. This product does not replace the need for caregiver supervision.



Periodic Safety and Maintenance checks:

To ensure user safety, regularly inspect the following components:

- Proper functioning of the buckles or hook and loop closure
- Signs of wear or damage to padding, stitching, or webbing
- Correct positioning of the product on the user

If any malfunctioning or wear is found on the belt, please contact your supplier for repair or replacement.

▲ In the event of a serious incident, report it to the manufacturer and the relevant local authority.

Cleaning instructions:

Hand wash or machine wash at 30°C. Do not bleach, tumble dry, or iron. To protect the product during washing, use a laundry bag. If required for more intensive cleaning, the product may be machine washed at 60°C.



Recycling:

At the end of its lifecycle, the product must be dismantled and disposed of in accordance with applicable local regulations for each material.

Warranty:

The product includes a 12-month warranty covering manufacturing defects.

Mounting a chest belt


The mounting of a 2-point hip belt should be determined based on the clinical goal. Always consult a seating specialist to determine optimal positioning.

A higher placement on the trunk provides increased control and stability, whereas a lower placement allows for greater freedom of movement.



Installation steps:

1. Identify the appropriate mounting points on the backrest of the wheelchair or the base of the seating system
2. Attach the mounting hardware securely to these points
3. Center the belt on the user
4. Thread the webbing through the mounting points, and adjust the length accordingly. Ensure the belt is tightened sufficiently to provide a secure and snug fit

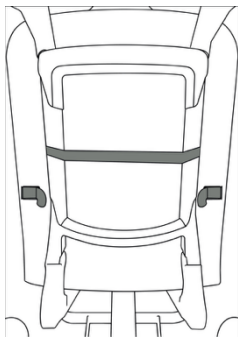

Want to know more about positioning?
Watch our e-learning.

E-learning





Adjustment of the belt

During use the belt can be adjusted using the ladder buckles. Ensure the belt remains centered on the user during adjustments.

Adjustment of the Chest Belt V1 is performed by tightening or loosening the hook-and-loop closure until the desired level of snugness is achieved.

Adjustments should only be made to:

- Facilitate easy fastening and unfastening
- Ensure the belt remains snug
- Fastening and loosening

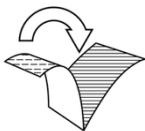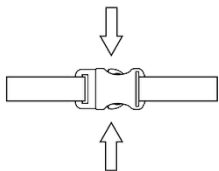

Symbols list



The product is manufactured by:
Kinetic Balance



Address:
Scheepmakerstraat 27
2222AB
Katwijk aan Zee
The Netherlands

Email: mail@positioningproducts.nl

Phone: +31 (0)71-2048424

Website: www.positioningproducts.nl

Chambre of commerce: 27339364

