



Safety & user manual

CHEST HARNESSES

Articles: 5.70001.X, 5.70002.X,
5.71032.X 5.71034.X,
5.81001.X, Z5.71022.X.

KBSUM-CH0925



This manual provides essential information for the safe use and maintenance of the chest harnesses range manufactured by Kinetic Balance. It is important that both the user and/or caregiver read and fully understand these instructions before use.

Intended use:

⚠ This product must be used in combination with a hip belt.

⚠ This product is exclusively designed for positioning the trunk of individuals seated in wheelchairs or other seating systems. It must **not** be used as a transportation belt or as a restraint device. Improper use may result in serious safety risks.

⚠ The belt should always be securely fitted around the trunk during use. Installation and adjustments must be carried out by a qualified seating specialist.

⚠ If a user's physical or cognitive abilities limit their ability to operate the belt safely, a caregiver must be present at all times. Caregivers should be properly instructed in fastening and releasing the product. This product does not replace the need for caregiver supervision.



Periodic Safety and Maintenance checks:

To ensure user safety, regularly inspect the following components:

- Proper functioning of the buckles
- Signs of wear or damage to padding, stitching, or webbing
- Correct positioning of the product on the user

If any malfunctioning or wear is found on the belt, please contact your supplier for repair or replacement.

▲ In the event of a serious incident, report it to the manufacturer and the relevant local authority.

Cleaning instructions:

Hand wash or machine wash at 30°C. Do not bleach, tumble dry, or iron. To protect the product during washing, use a laundry bag. If required for more intensive cleaning, the product may be machine washed at 60°C.



Recycling:

At the end of its lifecycle, the product must be dismantled and disposed of in accordance with applicable local regulations for each material.

Warranty:

The product includes a 12-month warranty covering manufacturing defects.

Mounting Kit:

SIZE	WEBBING CLIP			
	15mm	20mm	25mm	38mm
XXS	4pc			
XS		4pc		
S		4pc		
M			4pc	
L			4pc	
XL			4pc	

SIZE	CAM BUCKLE			
	15mm	20mm	25mm	38mm
XXS			4pc	
XS			4pc	
S			4pc	
M			4pc	
L			4pc	
XL			4pc	



SIZE	SLIDER BUCKLE			
	15mm	20mm	25mm	38mm
XXS	4pc			
XS		4pc		
S		4pc		
M			4pc	
L			4pc	
XL			4pc	



Mounting a Chest Harness


The mounting of a chest harness should be determined based on the clinical goal. Always consult a seating specialist to determine optimal positioning and belt angle.

Installation steps:

1. Determine the appropriate positioning of the chest harness on the user's trunk
2. Identify the correct mounting points for the upper straps, taking into account the user's shoulder width and seating position
3. Identify the mounting points for the lower straps so they align with the lower ribcage
4. Securely attach the mounting hardware to the identified points



5. Ensure that the pull direction of the upper straps runs straight toward the rear or at a maximum upward angle of 5°
6. Adjust and tighten the harness until it provides a secure and snug fit for the user

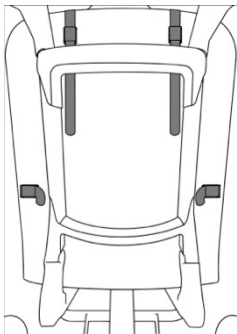

Want to know more about positioning?
Watch our e-learning.

E-learning




Adjustment of the harness

During use the harness can be adjusted using the D-ring(s). Ensure the harness remains centered on the user during adjustments.

Adjustments should only be made to:

- Facilitate easy fastening and unfastening
- Ensure the belt remains snug

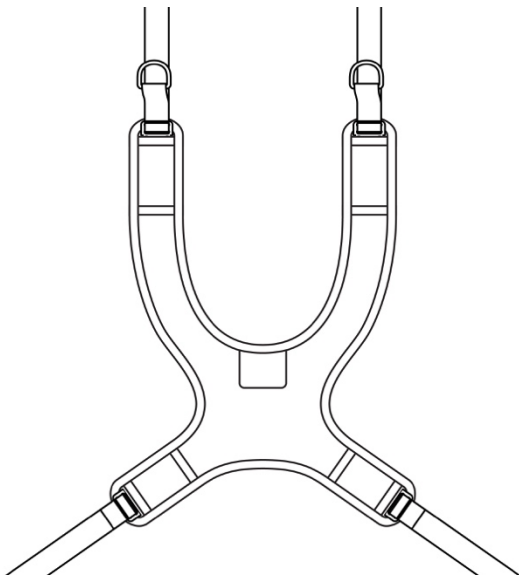


Flex control (Fitted Exclusively on non-stretch harnesses.)

5.71032.x, 5.71034.x, 5.81001.x

Use the Flex Control mechanism to alternate between the standard fit and a tighter configuration. Pull the strap toward the rear to achieve a tighter fit. To release the increased tension, pull the loop located at the top of the Flex Control.

Fastening and Loosening buckles

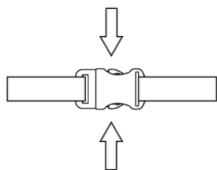

Symbols list



The product is manufactured by:
Kinetic Balance



Address:
Scheepmakerstraat 27
2222AB
Katwijk aan Zee
The Netherlands

Email: mail@positioningproducts.nl
Phone: +31 (0)71-2048424
Website: www.positioningproducts.nl
Chambre of commerce: 27339364

