



Safety & user manual

MULTIPURPOSE BAND

Articles: 500X

KBSUM-MPB0925



This manual provides essential information for the safe use and maintenance of the multipurpose band manufactured by Kinetic Balance. It is important that both the user and/or caregiver read and fully understand these instructions before use.

Intended use:

⚠ When the product is applied to the trunk, it must always be used in combination with a hip belt.

⚠ This product is exclusively designed for supporting the trunk and or lower or upper legs of individuals seated in wheelchairs or other seating systems. It must **not** be used as a transportation belt or as a restraint device. Improper use may result in serious safety risks.

⚠ The belt should always be securely fitted around the body during use. Installation and adjustments must be carried out by a qualified seating specialist.

⚠ If a user's physical or cognitive abilities limit their ability to operate the belt safely, a caregiver must be present at all times. Caregivers should be properly instructed in fastening and releasing the product. This product does not replace the need for caregiver supervision.



Periodic Safety and Maintenance checks:

To ensure user safety, regularly inspect the following components:

- Proper functioning of the hook and loop closure
- Signs of wear or damage to padding, stitching, or webbing
- Correct positioning of the product on the user

If any malfunctioning or wear is found on the belt, please contact your supplier for repair or replacement.

▲ In the event of a serious incident, report it to the manufacturer and the relevant local authority.

Cleaning instructions:

Hand wash or machine wash at 30°C. Do not bleach, tumble dry, or iron. To protect the product during washing, use a laundry bag. If required for more intensive cleaning, the product may be machine washed at 60°C.



Recycling:

At the end of its lifecycle, the product must be dismantled and disposed of in accordance with applicable local regulations for each material.

Warranty:


The product includes a 12-month warranty covering manufacturing defects.

Mounting a Multipurpose band

The mounting of a multipurpose band should be determined based on the clinical goal. Always consult a seating specialist to determine optimal positioning and band.

Installation steps:

1. Identify the appropriate location on the user
2. Tighten the band to ensure a snug and secure fit
3. use the hook-and-loop closure to secure the belt in place

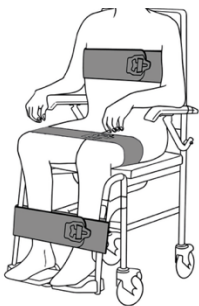

Want to know more about positioning?
Watch our e-learnings.

E-learnings



Adjustment of the belt

During use the belt can be adjusted using the hook-and-loop closure. Ensure the belt remains centered on the user during adjustments.



Symbols list



The product is manufactured by:
Kinetic Balance



Address:
Scheepmakerstraat 27
2222AB
Katwijk aan Zee
The Netherlands

Email: mail@positioningproducts.nl

Phone: +31 (0)71-2048424

Website: www.positioningproducts.nl

Chambre of commerce: 27339364

